


# Ana-Box™ Analog Sensor Interface-\$175

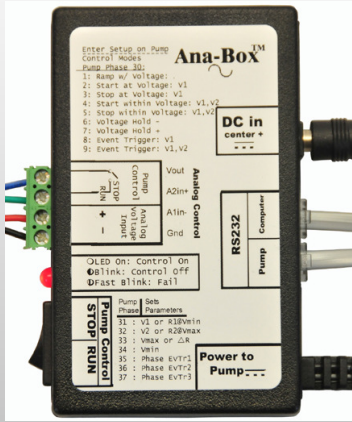

## MARKET EXCLUSIVE


Interface your pressure sensors, PLCs, sensors, joysticks, etc. with the Ana-Box™ Analog Sensor Interface to your pump.

*So capable...  
yet so incredibly simple to use!*


### Sample Ana-Box™ Operational

- Over 10 modes of operation
- Vary the pumping rate in proportion to a variable voltage input
- Start or stop pump at set voltages
- Vary the pumping rate to maintain a pressure or pH level  
I.e. maintain a pressure or pH level
- Attaches to any pump in the NE-1000 or NE-9000 Series
- Operates stand alone. Does not need a computer to operate
- Trigger program events according to voltage
- Optionally outputs event data to a computer for data logging
- No complicated user interface.
- OEM customizing available


**Ramping (gradient) with Voltage**  
Automatically adjusts the pumping rate proportionally to the voltage

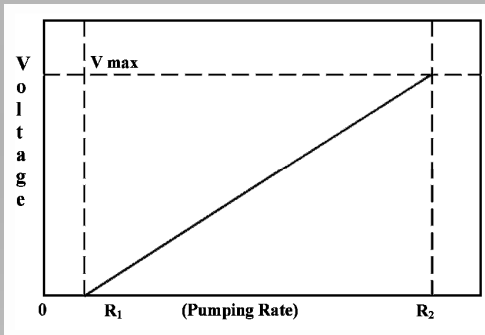


**Pressure Controlled Infusion**  
Automatically varies the pumping rate to maintain a pressure by utilizing your pressure sensor with a voltage output.

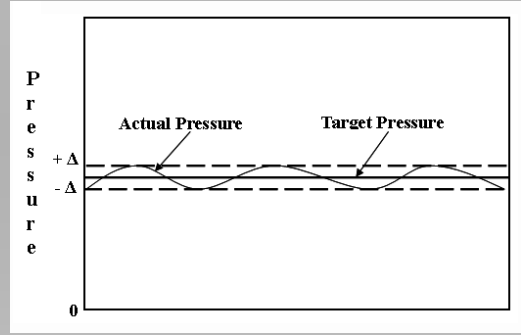


**Start Pump at Voltage Trigger**

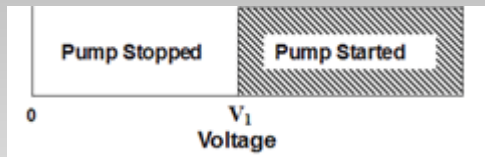

**Start Pump within Voltage**

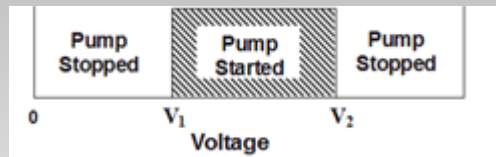

**Stop Pump at Voltage Trigger**

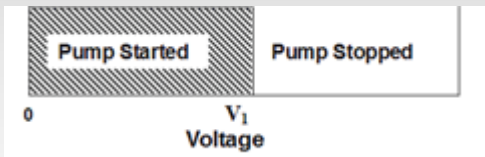

**Stop Pump within Voltage**

